
amieo

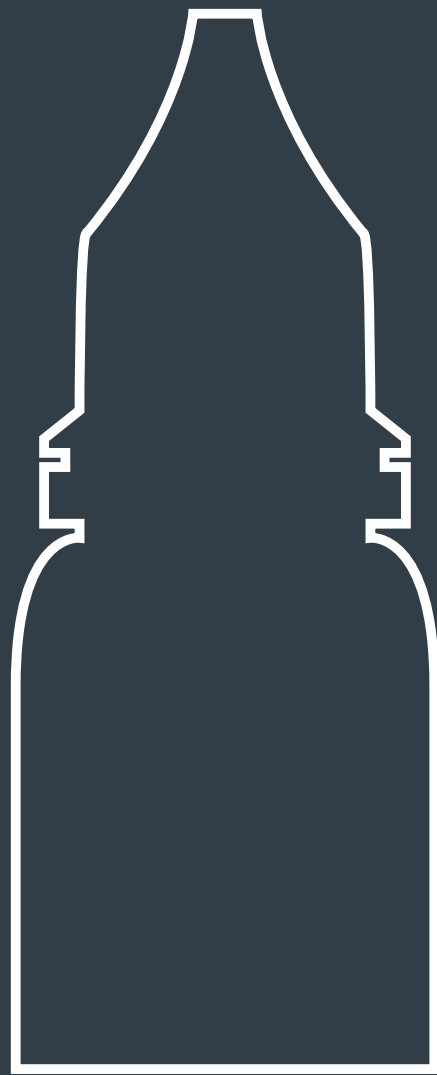MICROPIGMENTATION COLOR

DE

Inhaltsstoffe

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EN

Ingredients

The ingredients of each particular micropigmentation colors are listed together for all languages in the last section of this operating manual.

Contains nickel. Can cause allergic reactions.

Contains chromium (VI). Can cause allergic reactions.

Intended use

Mixture for use in tattoos or in permanent make-up. The micropigmentation color is implanted into the dermis with the aid of a needle cartridge and a micropigmentation device (permanent make-up).

The most common applications for micropigmentation are the following:

- Drawing fine hairs or shadowing of eyebrows
- Drawing an eyelid line
- Lip shading or lip contours
- Body tattoos

Contraindications **⚠ WARNING**

Micropigmentation is **not** permitted under the following circumstances:

- Hemophilia or other blood coagulation disorders
- Current consumption of blood thinners (e.g. acetylsalicylic acid, heparin, aspirin, warfarin)
- Uncontrolled diabetes mellitus
- Any form of active acne in the treatment area
- Dermatitis (e.g. skin tumors, keloids or extreme inclination for keloid formation, solar keratosis, warts and/or moles) in the treatment area
- Open wounds and/or eczema and/or skin rashes in the treatment area
- Scars in the treatment area
- Systemic infections and infective diseases (e.g. hepatitis A, B, C, D E or F; HIV infection) or acute local skin infections (e.g. herpes, rosacea)
- For eyelid pigmentation: Conjunctivitis
- During chemotherapy, radiotherapy or high-dosage corticosteroid therapy (recommendation: less than 4 weeks before the beginning until 4 weeks after the end of therapy)
- Up to 12 months after plastic surgery in the treatment area
- Up to 6 months after filler injections in the treatment area
- Allergy to topical anesthetics (local anesthesia)
- Under the influence of alcohol and/or drugs
- Pregnancy and lactation

Treatment of mucous membranes and eyeballs is strictly prohibited.

If there is the slightest doubt about customer safety, do not commence treatment or terminate treatment immediately.

The treatment must be terminated immediately in the event of:

- Excessive sensations of pain
- Fainting / dizziness
- Excessive release of wound fluid

As a general rule, recently treated areas of skin should be protected from UV radiation and sunlight.

In rare cases, the following problems may be triggered by the introduction of micropigmentation color:

- Color differences
- Pigment loss
- Allergic reactions to components of the aseptic pigment color

Use of micropigmentation colors **⚠ CAUTION**

It is critical to always observe the following notes prior to using the micropigmentation colors:

- Micropigmentation colors must only be used by medically qualified staff or staff qualified and trained for micropigmentation.
- The micropigmentation colors may only be used for the purposes listed under Intended Use.

-  **Shake the bottle very well before use.**
-  Hold the bottle diagonally at an angle of about 45° when removing the micropigmentation color.
- If the mixing balls are blocking the opening, firmly tap on the bottom of the bottle to loosen the mixing balls.
- Do not mix with other manufacturers' products.
- Perform a tolerance test before use *)
- Avoid contact with eyes. Rinse with water if the product comes in contact with the eyes.
- Inform the customer of potential risks.
- We strongly recommend the use of needle cartridges from MT.DERM, because hygienic results can only be guaranteed this way.
- Do not use the micropigmentation color if there is any reason to assume that the micropigmentation color is no longer in perfect condition (for example, as a result of transport damage or improper storage).
- Observe all national hygiene regulations.
- Make sure that any micropigmentation color remaining in the ink bottle is not contaminated during treatment.
- After removing the micropigmentation color from the ink bottle, do not fill it back in.
- Generally, store the ink bottle in a cool and dark location.
- Do not expose the ink bottle to direct sunlight over longer periods of time.

DE

EN

FR

IT

ES

PT

NL

RU

PL

CS

SK

EL

BG

RO

HR

SL

HU

SV

DA

FI

ET

LV

LT

MT

GA

- Only use a cloth soaked in disinfectant to clean the ink bottle.
- Do not drink or inhale the micropigmentation color.
- Use the micropigmentation color before the expiration date printed on the product label.
- Keep the micropigmentation color out of the reach of children.
- Remaining ink can be disposed of in household garbage.

* For the tolerance test, insert a small amount of the micropigmentation color to the sensitive skin behind the ear with a sterile needle. Wait 72 hours to ensure that there is no allergic reaction. The same test should then be performed on the skin area to be treated. Perform this test for each micropigmentation color to be used.

Risks **⚠ CAUTION**

- The micropigmentation color to be used remains in the skin and may in some cases not be able to be removed.
- Even though the color pigments it contains are very stable, the color may fade or change after use.
- There is no guarantee that the micropigmentation color will fade out entirely.
- Depending on the application by the trained staff and the condition of the skin, the color pigments may be deposited irregularly.
- Adverse allergic reactions may occur, even though no side effects have been reported so far.
- In some cases, it may not be possible to remove the micropigmentation color from clothing (or other materials).

Important notes **⚠ CAUTION**

Even though the product meets the most stringent requirements relating to quality and purity, the use of this product may trigger allergic reactions for some people due to their individual sensitivity to one or more ingredients in the color. The manufacturer will not assume any liability for allergic reactions of individual persons to colors, neither on the basis of an assurance of serviceability or expediency nor for any other reason. If the tolerance of the color is unknown from former applications, perform a color tolerance test prior to use or consult with a dermatologist.

⚠ WARNING

The manufacturer assumes no liability for the consequences of attempts to remove the colors. In particular, the manufacturer will not assume any liability for complications resulting from the removal of micropigmentation color by way of any thermal, mechanical or other process (such as laser treatment).

For any complaints, generally specify to the lot code printed on the product label.

Hazard warnings

Flammable liquids Category 3, Warning!
H226 „Flammable liquid and vapor“.

Definition of the symbols used on the label and product packaging

The following symbols can be used on the product label and the packaging for identification purposes:

Manufacturer

See instruction for use for further information

Catalog number, order number

Lot code

Effective until at least

Best used within twelve months after opening

DE

EN

FR

IT

ES

PT

NL

RU

PL

CS

SK

EL

BG

RO

HR

SL

HU

SV

DA

FI

ET

LV

LT

MT

GA

Manufacturer's Quality Declaration

- MT.DERM GmbH is ISO 13485 certified.
- The ingredients of the micropigmentation colors meet the requirements of Annex XVII of the (EC) Regulation no. 1907/2006 (REACH) concerning substances in tattooing colors or permanent make-up.
- The microbiological harmlessness of all micropigmentation colors has been confirmed by an accredited laboratory.
- All micropigmentation colors are free of preservatives.
- All micropigmentation colors are backed by a seal of authenticity.

Warning symbols and statements used in this Instructions for use

WARNING

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury

CAUTION

CAUTION indicates a hazardous situation which, if not avoided, could result in minor or moderate injury

Manufacturer

MT.DERM GmbH

Blohmstraße 37-61

12307 Berlin, Germany

Tel: +49 30 7676 6220-0

www.mtderm.de, info@mtderm.de

DEEP, EVO2, Evolution Line, Organic Line and Scalp Micropigmentation colors - package insert

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GA

Declaration of Conformity (DoC)

REACH Annex XVII according to Commission Regulation (EU) 2020/2081

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declares under sole responsibility that the mixtures listed in the above section for use in tattoos or permanent make-up comply with the provisions of the following regulation:

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COMMISSION REGULATION (EU) 2020/2081 of 14 December 2020 amending Annex XVII to Regulation (EC) No. 1907/2006 of the European Parliament and of the Council concerning the Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) as regards substances in tattoo inks or permanent make-up.

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This Statement is made for the manufacturer by:

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© 06.2024
The manufacturer: MT.DERM GmbH